



Always on hand  
to help with hygiene –  
**in Europe and around  
the world.**



neodisher® –  
for maximum  
**Hygiene  
Safety.**

Chemische Fabrik  
Dr. Weigert GmbH & Co. KG  
Mühlenhagen 85  
20539 Hamburg  
Tel.: +49-40-7 89 60-0  
info@drweigert.de

Dr. Weigert  
Handelsgesellschaft m.b.H.  
Wienerbergstraße 11/12 a  
1100 Wien  
Tel.: +43-1-99460 6275  
info@drweigert.at

Dr. Weigert (Schweiz) AG  
Baarerstrasse 14  
6300 Zug  
Tel.: +41-41 561 32 90  
info@drweigert.ch

Dr. Weigert Nederland BV  
Narcisstraat 14  
9404 RK Assen  
Tel.: +31-592-31 93 93  
info@drweigert.nl

Dr. Weigert UK Ltd  
Watling Court, Orbital Plaza,  
Watling Street, Cannock  
WS11 0EL  
Tel.: +44-1543 404 633  
enquiry@drweigert.com

Dr. Weigert Polska Sp. z o.o.  
ul. Wybrzeże Gdynskie 6A  
01-531 Warszawa  
Tel.: +48-22-616 02 23/31  
office\_pl@drweigert.com

Dr. Weigert España  
Edificio Cuzco IV  
Paseo de la Castellana 141,  
Planta 8  
28046 Madrid  
Tel.: +34-91 572 65 77  
info.spain@drweigert.com

Dr. Weigert France SAS  
22 avenue des Nations  
Immeuble le Raspail  
BP 88035 Villepinte  
95932 Roissy CDG Cedex  
Tel.: +33-1-48 67 90 33  
info@drweigert.fr

Dr. Weigert Hungária Kft.  
Hunyadi János út 16  
1117 Budapest  
Tel.: +36-1-237 06 04  
info@drweigert.hu

GB 1023 / 11-2013 © Dr. Weigert, Hamburg

Dr. Weigert management system – Certified according to DIN EN ISO 9001/13485/14001

[www.drweigert.com](http://www.drweigert.com)


Reprocessing of Bedpans  
and Urine Bottles.

# neodisher® SBN plus – Special Rinse Aid for Bedpans.

## neodisher® SBN plus

Effective special rinse aid for the reprocessing of bedpans and urine bottles

- With fresh linden blossom fragrance against unpleasant odours
- For washer disinfectors with steam disinfection

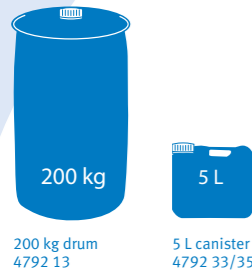

\*The product is available in two different 5 L canisters; the canister with art. no. 4792 35 is higher, the canister with art. no. 4792 33 is flatter and oblong.

### Fields of application:

- For water softening and rinsing when reprocessing bedpans and urine bottles in bedpan washer disinfectors with steam disinfection
- For fast and stain-free drying
- Suitable for acid-compatible bedpan washers
- For removing stubborn soiling we recommend the additional use of **neodisher® SBR extra**, a special detergent for bedpans and urine bottles

### Disinfection with hot steam:

- With bedpan washers using a thermal reprocessing process the disinfecting step is effected with a hot water-steam mixture
- The steam generator is the heart of the bedpan washer and must be protected against the formation of lime which may cause longer heating times, higher power consumption, overheating and thus could lead to irreparable damage to the heating elements



## Our Competence – Your Benefit:

- acidic special rinse aid based on organic acids
- with special softening components for effectively complexing water hardness
- with pleasant fresh linden blossom fragrance
- with surfactants in a low-foam formula



- effectively prevents the formation of lime in the steam generator, in the flushing chamber and on the item to be reprocessed
- suitable for all degrees of water hardness
- reliably covers unpleasant odours
- for optimized, rapid and stain-free drying without lime deposits



Pay attention to acid compatibility of dosing systems and steam generators. For further information on products and applications, please refer to the respective product data sheet.

# neodisher® Cleaning – Water Softening & Rinsing.

## neodisher® SBR extra

Material-compatible detergent, even removes stubborn residues

- Improves the cleaning performance, removes particularly stubborn residues (such as lime soap stains and residues from cytostatics in faeces)
- Suitable for bedpan washers equipped with a respective dosing pump




200 kg drum 4007 13  
5 L canister 4007 33/35\*

## neodisher® SBK

Combined softening and rinsing agent based on organic acids

- Prevents the formation of lime in the steam generator, pipes and jets
- For optimum stain-free drying
- Suitable for acid-compatible bedpan washers

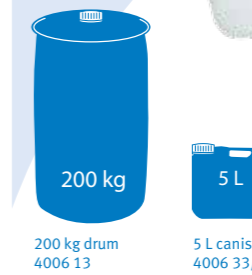


200 kg drum 4006 13  
5 L canister 4006 33/35\*

## neodisher® SBN

Combined softening and rinsing agent based on organic acids and surfactants

- Prevents the formation of lime in the steam generator, pipes and jets even when using hard water
- For optimum stain-free drying
- Suitable for acid-compatible bedpan washers

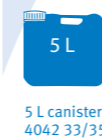

5 L canister 4042 33/35\*

## neodisher® SBX

Mildly-alkaline softening agent, complexes the water hardness

- Prevents the formation of lime deposits in the steam generator
- Suitable for all bedpan washers



5 L canister 4015 33/35\*

\*The product is available in two different 5 L canisters; the canister with art. no. xxxx 35 is higher, the canister with art. no. xxxx 33 is flatter and oblong.